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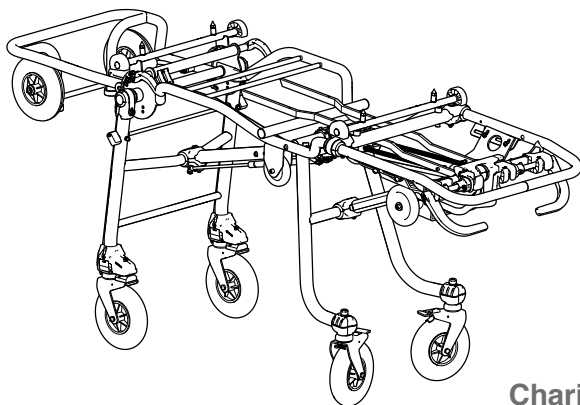
EN

1.	MODELS	15
2.	INTENDED USE	15
3.	REFERENCE STANDARDS	15
4.	INTRODUCTION	15
4.1	Using the instructions for use	15
4.2	Device labelling and traceability	16
4.3	Symbols	16
4.4	Warranty and service	16
5.	WARNINGS	16
6.	SPECIFIC WARNINGS	17
6.1	Physical requirements for operators	17
7.	RESIDUAL RISK	18
8.	TECHNICAL DATA AND COMPONENTS	18
9.	COMMISSIONING	19
10.	OPERATING CHARACTERISTICS	19
11.	PROPER USE	20
11.1	Transport vehicle requirements	20
11.2	Parking brakes	20
11.3	Twist system activation	20
11.4	Lowering the stretcher	20
11.5	Lifting the stretcher	21
11.6	Loading and unloading the stretcher from the vehicle	21
11.7	Adjusting the stretcher length	22
11.8	Method of use for coffin support	22
12.	CLEANING AND MAINTENANCE	23
12.1	Cleaning	23
12.2	Ordinary maintenance	23
12.3	Periodic overhaul	23
12.4	Special maintenance	24
12.5	Life span	24
13.	TROUBLESHOOTING TABLE	24
14.	ACCESSORIES	24
15.	SPARE PARTS	24
16.	DISPOSAL	24

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1. MODELS

The following basic models may be subject to implementation or change without notice.

- ZCROSS – ROLL-IN STRETCHER
- ZCROSS – ROLL-IN STRETCHER WITH TRAY
- ZCROSSOVER – ROLL-IN STRETCHER
- ZCROSSOVER – LOW ROLL-IN STRETCHER
- ZCROSSOVER – ROLL-IN STRETCHER WITH TRAY
- ZCROSSOVER – LOW ROLL-IN STRETCHER WITH TRAY

The models are available in standard and low versions.

The above-listed models may be marketed in different variants and, depending on the model, may include accessories.

2. INTENDED USE

Spencer Z-Cross and Z-Crossover stretchers are compact and lightweight roll-in stretchers with variable heights for the transport of bodies. The legs can be folded independently for loading by engaging the corresponding controls and they open automatically during unloading. The device has a safety system on both the front and rear legs that is automatically activated when the stretcher is loaded onto the transport vehicle. Large diameter wheels, which are fixed at the front and self-aligning at the rear, improve loading and transport.

USERS AND INSTALLERS

Z-Cross and Z-Crossover stretchers are devices intended for professional use only. All operators must be trained to transport bodies safely and efficiently. For body loading techniques for especially heavy subjects, for operations on steep terrain or in special and unusual circumstances, the presence of multiple operators is recommended (not just 2 as expected in standard conditions).

CONTRAINDICATIONS AND UNWANTED SIDE EFFECTS

No particular contraindications or side effects are known with relation to use of the device, as long as it is used in accordance with the instructions for use.

USER TRAINING

Note: despite all efforts, laboratory tests, trials, and instructions for use, standards do not always reproduce practice, so the results obtained under actual conditions of product use in the natural environment may sometimes differ significantly.

The best instructions are the continuous practice of use under the supervision of competent and trained personnel.

- Regardless of your level of experience with similar devices in the past, you should carefully read and understand the contents of this manual before installing, operating, or servicing this product. In case of any questions, please contact Spencer Italia S.r.l. for the necessary clarifications.
- The product must be used only by personnel trained in the use of this product and not on other similar products.
- The suitability of the users for use of this product can be attested by the training registration, in which trained persons, trainers, date and place are specified. **This documentation must be kept for at least 10 years after the end of the product's life and must be made available to the competent authorities and/or the Manufacturer when requested. In the absence of such documentation, the relevant bodies will apply any foreseen sanctions.**
- Do not allow untrained persons to help while using the product, as they may cause injury to themselves or others.
- The product must be put into use only by personnel trained in the use of this product and not on other similar products.

Note : Spencer Italia S.r.l. is always available for training courses.

INSTALLER TRAINING

Device installation must be carried out by qualified personnel who have been trained and certified by Spencer Italia S.r.l.

The installer must strictly adhere to these instructions as well as to the state of the art for vehicle installations.

3. REFERENCE STANDARDS

EU Regulation 2023/988

General product safety regulation

4. INTRODUCTION

4.1 USING THE INSTRUCTIONS FOR USE

The purpose of these instructions is to provide professionals with the information necessary for safe and appropriate use and maintenance of the device.

Note: The Manual is an integral part of the device and therefore it must be kept for the entire life of the device and must accompany it in any changes of use or ownership. If any instructions for use for products other than the one received are present, please contact the Manufacturer immediately before use.

The instructions for use of Spencer products can be viewed and downloaded from the Spencer website at www.spencer.it/en/emergency-medical-solutions/user-manuals or by contacting the Manufacturer. Exceptions are items whose essential nature and reasonable and foreseeable use are such that instructions are not necessary, in addition to the following warnings and instructions on the label.

Regardless of your level of experience with similar devices in the past, it is advisable to carefully read and understand the contents of this manual before installing, operating, or servicing this product.

IT

EN

DE

FR

ES

PT

BG


HU

LT

EL

4.2 DEVICE LABELLING AND TRACEABILITY

Each device is provided with a label, placed on the device itself and/or on the packaging, which contains the Manufacturer's identification data, product, serial number (SN) or lot number (LOT). This must never be removed or covered.

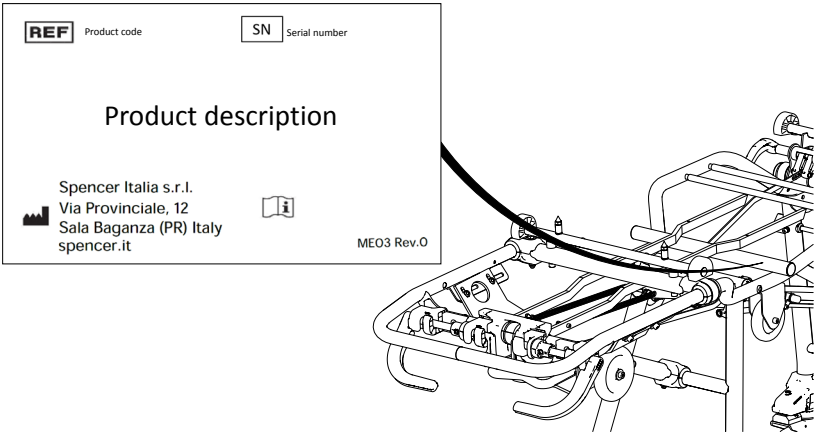In the event of damage or removal, request a duplicate from the Manufacturer, or else the warranty will be void as the device can no longer be traced.

4.3 SYMBOLS

Symbol	Meaning	Symbol	Meaning
	Manufacturer		Danger – Indicates a hazardous situation that may result in a situation directly related to serious injury or death.
	Date of manufacture		See the user manual
	Product code		Serial Number

4.4 WARRANTY AND SERVICE

Spencer Italia S.r.l. guarantees that products are free from defects for a period of **one year from the date of purchase**.
For information regarding correct interpretation of the instructions for use, maintenance, installation or return, please contact Spencer After Sales tel. +39 0521 541111, fax +39 0521 541222, e-mail service@spencer.it.

To facilitate service, always indicate the lot number (LOT) or serial number (SN) on the label attached to the package or device itself.

Warranty and service conditions are available at www.spencer.it/en/support/terms-and-conditions

To ensure product traceability and protect the maintenance and servicing procedures for your devices, Spencer has made the following portal available to you www.spencer.it/en/services/report-an-issue

Note: Record and keep with these instructions: lot (LOT) or serial number (SN), if present, place and date of purchase, date of first use, date of checks, user name and comments.

5. WARNINGS

Warnings, notes, and other important safety information are provided in this section and are clearly visible throughout the manual.

At least every 6 months, it is important to check for updated instructions and any changes involving your product. This information is freely available on the website www.spencer.it on the specific product page.

PRODUCT FEATURES

Use of the product for any purpose other than that described in the User Manual is prohibited.

- Before each use, always check the conditions of the product, as specified in the User Manual. In the event of faults/damage that could compromise its functionality/safety, immediately remove it from service and contact the Manufacturer.
- If the product is found to be malfunctioning, immediately use a similar device to ensure continuity of ongoing operations.
- The product must not be tampered with (modification, tweaking, additions, repair). Should these operations be performed, we decline any responsibility for incorrect operation or any damage caused by the product itself; furthermore, the product warranty will be void.
- Anyone who modifies or causes products manufactured by Spencer Italia S.r.l. to be modified or re-processed in such a way that they no longer serve the intended purpose or no longer provide the intended performance shall comply with the conditions valid for the first placing on the market.
- When using the devices, position and adjust them in such a way that they do not hinder operator works or the use of any other equipment.
- Be sure to take every precaution to avoid hazards from contact with blood or body secretions, if applicable.
- Always comply with the maximum capacity indicated in the User Manual. Maximum load capacity means the total weight distributed according to human anatomy. When determining the total weight load on the product, the operator should consider the weight of the body, equipment and accessories. In addition, the operator should assess whether the overall size of the body reduces the functionality of the product.
- Make sure that operators are in good physical condition before lifting, as listed in the User Manual.
- **The maximum weight, which weighs on each operator, must comply with local health and safety requirements.**
- Avoid contact with sharp objects.
- Device installation must be carried out by qualified personnel trained and enabled by Spencer Italia S.r.l. The times and methods of such shall be agreed upon between the customer and our Sales Offices.
- Operating temperature: from -10°C to + 50°C.

STORAGE

- The product must not be exposed or come into contact with thermal sources of combustion or flammable agents, but must instead be stored in a dry, cool place, away from light and sun.
- Do not store the product under other more or less heavy materials that may damage the structure of the product.
- Store and transport the product with its original packaging, otherwise the warranty shall be invalidated.
- Storage temperature: -20°C to +60°C.

REGULATORY REQUIREMENTS

As Distributor or End-User of the products manufactured and/or marketed by Spencer Italia S.r.l., users are strictly required to be familiar with the legal provisions in force in the country of destination of the goods, applicable to the devices to be supplied (including regulations relating to technical specifications and/or safety requirements) and, therefore, to understand the requirements necessary to ensure compliance of the products themselves with all legal requirements of the territory.

- Promptly and in detail notify Spencer Italia S.r.l. (already in the quotation request phase) about possible fulfillments by the Manufacturer necessary for the compliance of products with specific legal requirements of the territory (including those deriving from regulations and/or regulatory provisions of another nature).
- Act with due care and diligence to help ensure compliance with the general safety requirements of the devices placed on the market, providing end-users with all the information necessary to carry out periodic revisions on the supplied devices, exactly as indicated in the User Manual.
- **Participate in safety checks on products** placed on the market, transmitting information regarding product risks to the Manufacturer as well as to the Competent Authorities for their respective actions.
- Without prejudice to the above, the Distributor or End-User shall assume wider liability related to non-compliance with non-fulfilment of the above-mentioned obligations, with consequent obligation to indemnify and/or hold Spencer Italia S.r.l. harmless from any possible injurious effect.

6. SPECIFIC WARNINGS

To use the product, you must also have read, understood and carefully follow all the instructions in the user manual.

- Do not use if the device or parts of it are punctured, torn, frayed, or excessively worn.
- Before any type of handling, make sure that operators have a secure grip on the device.
- Do not drag the device on uneven surfaces.
- Do not lift with a crane or other mechanical lift.
- Do not use drying machines.
- The device is a transport device and **cannot be used as a stationing device**.
- Do not use with devices other than those expressly approved by the manufacturer.
- For body loading techniques for especially heavy bodies, for operations on steep terrain or in special and unusual circumstances, the presence of multiple operators is recommended (not just 2 as expected in standard conditions).
- Make sure that the body is adequately immobilized before loading it on the stretcher. A failure to immobilize can cause serious damage to the patient.
- Make sure that the sheet does not interfere with any stretcher handling or control mechanisms.
- Do not move the device if the weight is not properly distributed.
- Use only the perimeter frame to move the stretcher and not other points not intended for this purpose.
- Avoid excessive force when loading the stretcher into the vehicle: unnecessary force can cause damage and can adversely affect the fastening system.
- **Parking brakes are aids for the operator but they do not in any way replace his or her supervision.**
- Pay the utmost attention to any obstacles (water, ice, debris, etc.) present on the route, as they may cause the operator to lose balance and compromise proper functioning of the device. If you cannot clear the way, choose an alternative route.
- Condensation, water, ice and accumulations of dust can affect correct operation of the device, making it unpredictable and causing a sudden change in the weight the operators have to support.
- For height differences greater than 10mm, the device must be lifted, taking care to grasp it from the structure and not from other points not intended for this purpose.
- Once the wheels of the loading trolley are positioned on the ambulance platform, the wheels of the front leg must be between 5 cm and 10 cm from the ground, so that the front leg can be opened and locked securely. Check the height of the ambulance loading platform after each use; if altered, the trolley must be adjusted immediately by the manufacturer or a specialized technician authorized by the manufacturer. Failure to do so will result in the manufacturer declining all responsibility for the correct functioning of the device or for any damage caused by it.
- If the vehicle is equipped with pneumatic or hydraulic suspensions, the loading height must be adjusted taking into account the worst-case conditions and/or the working conditions foreseen by the outfitter.
- Any problems regarding the use and/or safety risks associated with this system are not the responsibility of the manufacturer.
- **Inadequate installation of the loading floor can lead to buckling and subsequent damage to the welds of the front legs.**
- **Improper installation of the loading floor can lead to abnormal operation of the device and cause harm to the user.**
- Do not alter or modify the stretcher to adapt it to the emergency vehicle: doing so could result in unpredictable operation and damage to the rescuers and shall void the manufacturer's warranty and release the Manufacturer from all liability.
- Do not operate the variable height system without carefully assessing the weight of the body with the patient and any accessories. Operators should be able to fully support the load when changing heights. Incorrect judgement may result in the stretcher falling unexpectedly and posing a risk to the operators.
- Do not place any magnetic parts between the stretcher and the fixing system as they may interfere with the stretcher's attachment and release systems.
- Do not place limbs and/or objects between the legs and the frame, in the vicinity of the leg lifting pistons or in general between moving parts as this may result in crushing injuries.
- If the twist system has been activated, make sure that the wheels are locked again before loading the stretcher into the vehicle as this can be risky and difficult with the front wheels unlocked.
- Movement with four swivel wheels can be very difficult on sloping or uneven ground. Carefully assess the conditions of use before unlocking the front wheel pivot.
- Use by untrained personnel may result in injury to the operator or third parties.
- Inadequate disinfection procedures may lead to risk of cross-infection.
- Partial opening of the legs may result in the device falling to the ground. Make sure that the legs are properly locked before any handling and that the pistons are fully extended and stable.
- If the stretcher is not locked on the fixing system or if it is not correctly positioned, this can result in dangerous movements, especially in the event of strong vehicle deceleration, which can cause damage to the operators. Always make sure that the locking system is properly engaged.
- Failure to observe the warnings for operators can lead to risks of crushing caused by the moving mechanisms.
- Accidental operation of the variable height system may result in the stretcher falling causing injury to the operators. Make sure that the unlocking handle is not accidentally activated.
- The operators should prepare to fully load the stretcher with the body and accessories before operating the variable height unlocking control. Operating this control without applying sufficient support force will result in the stretcher falling unexpectedly and causing injury to the operators.
- **Failure to read and understand the product instructions may result in injury to the operators.**

6.1 PHYSICAL REQUIREMENTS FOR OPERATORS

The device is intended for professional use only. Do not allow untrained persons to help while using the product, as they may cause injury to themselves or others.

Operators using the device should be physically able to use the device and have good muscle coordination, as well as strong back, arms, and legs for lifting and holding and be able to grasp the device firmly with both hands.

IT

EN

DE

FR

ES

PT

BG

HU

LT

EL

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Users must be able to safely lift and handle the weight of the stretcher and body together and any other equipment used with the device.
For body loading techniques for especially heavy bodies, for operations on steep terrain or in special and unusual circumstances, the presence of multiple operators is recommended (not just 2 as expected in standard conditions).
The ability of each operator must be assessed before the definition of rescuer roles in use of the stretcher.

7. RESIDUAL RISK

No residual risks have been identified, i.e., risks that could arise despite compliance with all the warnings contained in these instructions for use.


8. TECHNICAL DATA AND COMPONENTS

Note: Spencer Italia S.r.l. reserves the right to make changes to specifications without notice.


Definition of height under wheel depending on model:

Model	Height under wheel [cm]	
	Min	Max
Standard Cross / Crossover	58	67
Low ZCross / ZCrossover	47	55

Roll-in stretcher:

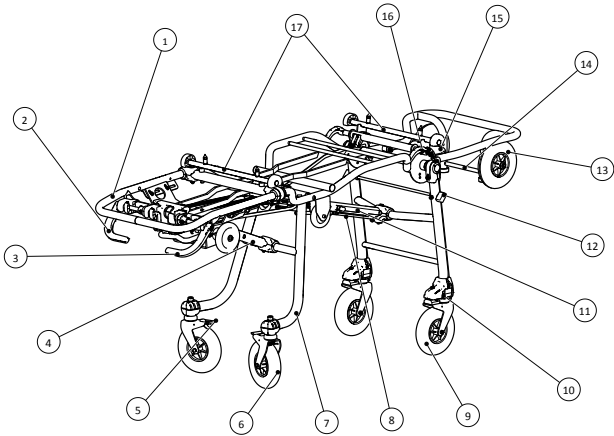Roll-in stretcher with tray:

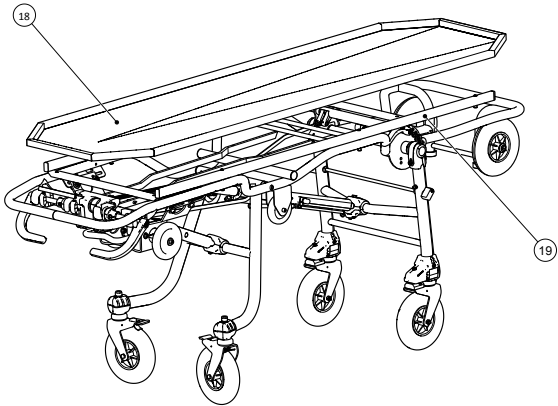No.	DESCRIPTION	MATERIAL
1	Perimeter frame	Steel
2	Rear leg unlocking handle (green)	Nylon
3	Front leg unlocking handle	Nylon

4	Rear piston	Steel
5	Rear wheel brake	Iron
6	Rear wheel with bracket	PU
7	Rear leg	Steel
8	Gas spring	Steel
9	Front wheel	PU
10	Twist system	Steel/Nylon
11	Front piston	Steel
12	Twist control lever	Nylon
13	Loading trolley wheels	PU
14	Crossbar for front fixing	Steel
15	Loading trolley unlocking control	PU/Steel
16	Variable height unlocking lever	Nylon
17	Mortuary trolley support crossbars	Steel/PA6
18	Tray support frame	Stainless steel
19	Tray	Stainless steel

FEATURE	DIMENSION
Length [mm]	1970
Width [mm]	560
Wheel diameter [mm]	Ø 200 ± 5%
Load capacity [kg]	170 kg (ZCross), 250kg (ZCrossover)
Weight without standard accessories [kg] (2)	38 (ZCross), 41 (ZCrossover), 48 (ZCross with tray), 51 (ZCrossover with tray)

(2) The weight may vary by no more than ±1 kg depending on the model.

9. COMMISSIONING

For first use, check that:

- Packaging is intact and has protected the device during transportation.
- Check that all parts included in the packing list are present
- General functionality of the device
- The vehicle is equipped with a Spencer fixing system dedicated to the stretcher
- The stretcher floor is properly levelled
- The stretcher floor is wide and long enough to accommodate the stretcher and its accessories without becoming an obstacle
- During loading and unloading from the vehicle, the front leg wheel must be between 5 cm and 10 cm from the ground to allow the front leg to open and lock safely - see image in point 11.7.
- The fastenings must keep the stretcher attached to the structure of the vehicle

Do not modify the stretcher structural, leverage or traction parts in any way as this could cause damage to the operators.

⚠ Failure to take the above measures will preclude safe use of the device, resulting in risk of damage to the operators and the device itself.


To facilitate insertion of the stretcher into the vehicle, eliminate sharp edges on the edge of the loading bed. The stretcher should be secured using Spencer hooks in such a way as to prevent any movement during transport in the vehicle, even in difficult driving conditions. Practice with a stretcher without a body before putting the device into service.

For subsequent use, perform the operations specified in paragraph 12.

If the above conditions are met, the device may be considered ready for use; otherwise, you must immediately remove the device from service and contact the Manufacturer.

Do not alter or modify the device arbitrarily, as doing so could result in unpredictable operation and harm to the operators and will void the warranty and release the Manufacturer from all liability.

10. OPERATING CHARACTERISTICS

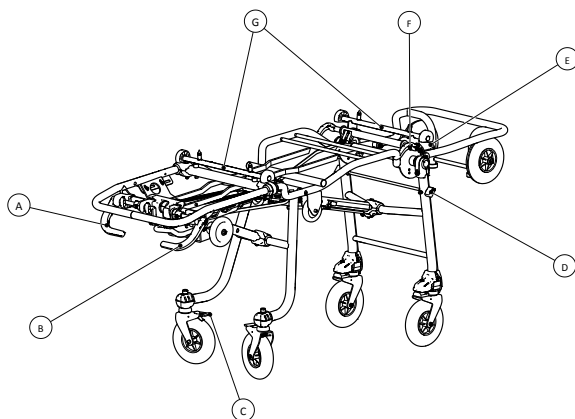IT

EN

DE

FR

ES

PT

BG

HU

LT

EL

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ELEMENT	DESCRIPTION	FUNCTION
A	Rear leg unlocking handle (green)	When engaged, it unlocks the movement of the front leg, allowing it to be closed.
B	Front leg unlocking handle (red)	When engaged, it unlocks the movement of the front leg, allowing it to be folded.
C	Rear wheel brake	When engaged, it brakes the rear wheels.
D	Twist control lever	When engaged, it also pivots the front wheels. This function may not be available on your model.
E	Variable height unlocking lever	Unlocks movement of the front leg in order to place the stretcher at intermediate height.
F	Loading trolley unlocking control	Unlocks the loading trolley and allows it to be raised in order to reduce the length of the stretcher when it has to be moved in confined spaces.
G	Support crossbars	Used for handling systems for transporting bodies.

11. PROPER USE

11.1 TRANSPORT VEHICLE REQUIREMENTS

The stretcher is designed to be loaded and unloaded from the compartment of a vehicle used for transporting bodies. The vehicle requirements should be:

- Level stretcher floor
- Stretcher floor that is wide and long enough to accommodate the stretcher without becoming an obstacle


The wheels of the front legs must be at a safe distance of at least 5 cm from the ground during loading/unloading or when the loading trolley is still on the vehicle, to allow the front leg to open safely.

 Failure to take the above measure will preclude safe use of the device, resulting in risk of damage to the operators and the device itself.

11.2 PARKING BRAKES

To engage the parking brakes, simply press on the strips on the rear wheel brackets with your foot. To disengage the parking brakes, simply press down on the opposite side of the parking brake pedal and they will return to their original position with a slight click.

Never leave the device unattended, even if the parking brakes have been engaged.

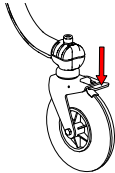
11.3 TWIST SYSTEM ACTIVATION

This system, if included in the model, allows the front wheels to pivot. To unlock rotation, turn the lever located on the side of the front leg counter-clockwise.


When rotation is no longer necessary, reposition the control in its original position by turning it clockwise and push the stretcher forward firmly; the wheels will then be locked again.

 When using the twist system, at least one additional operator must be present to steer the head end of the stretcher, as the stretcher is much more difficult to manoeuvre with all four wheels unlocked.

 Be sure to carefully assess the ground conditions before operating the twist system, as the stretcher may move in an unexpected and uncontrolled manner.

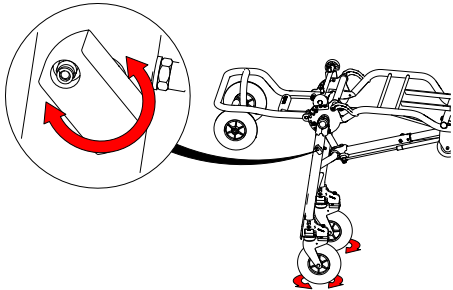11.4 LOWERING THE STRETCHER

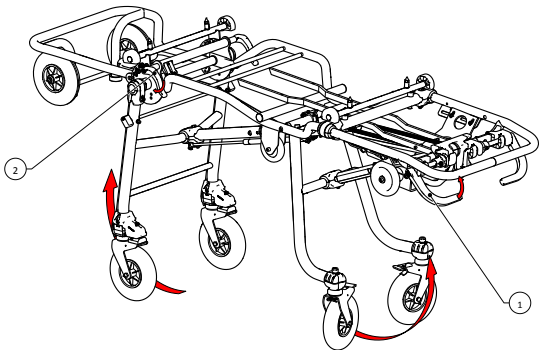To facilitate body transfer, it is advisable to move the stretcher to the intermediate height position or to the lowered position if there is no body present.

Procedure for loading/unloading a body onto the stretcher at intermediate height:

- If included in the model, engage the rear leg unlocking lever, hold it in place and lift the stretcher slightly to unlock the rear leg moving mechanism (no. 1 in figure). Lower

the stretcher about 10 cm, release the unlocking lever and support the stretcher until the intermediate height position is reached. Make sure that the stretcher has reached a stable position. Engage the brakes on the rear wheels.

- To lower the front part of the stretcher, engage the red lever (no.2 in figure), pushing the front leg slightly towards the loading trolley. The weight of the stretcher, body and any equipment attached to the stretcher must be supported. After starting the lowering motion, release the handle while keeping your grip on the frame until the stretcher reaches the intermediate height position. Make sure that the position reached is stable.


Lowering the stretcher to the ground

This procedure allows the stretcher to reach the minimum distance from the ground.

⚠ This manoeuvre cannot be carried out when a body is on the device.

- Raise the stretcher from the side of the feet until the loading trolley wheels rest on the ground.
- With the loading trolley wheels resting on the ground and ready to bear the weight of the stretcher, engage both leg unlocking levers, push towards the stretcher to unlock the front leg moving system and then the stretcher towards the ground. The stretcher is now in the lowest position.

Note: The parking brakes will not perform their function when the stretcher is fully lowered. Make sure that the stretcher is held in place by at least one operator.

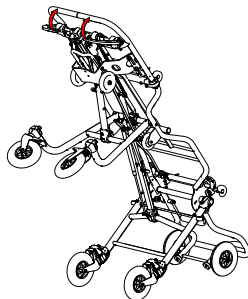
11.5 LIFTING THE STRETCHER

To bring back the stretcher back to the standard height, starting from any of the previous configurations, the operators must coordinate while simultaneously lifting the front and rear of the stretcher. Only carry out these operations after making sure that the body is well immobilised.

- At the rear, grasp the part of the frame near the feet of the body.
- At the front, similarly, grab the frame above the loading trolley or, for stretchers with detachable stretchers, grab the telescopic handles after making sure the stretcher is properly secured to the trolley.
- Raise the assembly until the locking mechanisms are properly engaged.


⚠ Only use the frame for lifting.

- Do not use the boards or other areas that are not designed for this purpose.

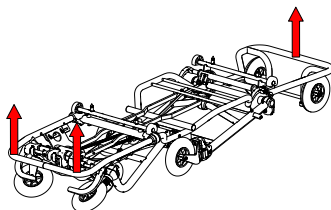11.6 LOADING AND UNLOADING THE STRETCHER FROM THE VEHICLE

⚠ Stretcher loading and unloading are among the most delicate steps during the use of stretchers. It is necessary to follow all the warnings in this manual very carefully, practising under simulated conditions before putting the product into service. The body must always be properly immobilized. Always ensure that the conditions of use allow the front legs to open properly before any loading and/or unloading operations.

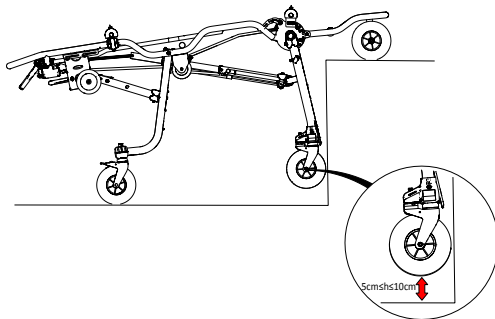TO LOAD THE STRETCHER ONTO THE VEHICLE:

- The operator performing the loading should align with the stretcher and the fixing system so that both devices are in front of him/her and in the intended alignment for fixing.
- Move the wheels of the loading trolley forward into the vehicle until the front legs of the stretcher come into contact with the bodywork of the vehicle.
- In the event of a slide/ramp, the wheel of the loading trolley should be in a horizontal position before the front legs come into contact with the vehicle bodywork.
- During loading/unloading, i.e. when the loading trolley is still resting on the vehicle, the wheels of the front legs must be at a safe distance of between 5 cm and 10 cm from the ground to allow the front leg to open safely.
- Engage the right red lever (letter B par.10) to unlock movement of the front legs while continuing to push the stretcher inside the vehicle until the rear legs are brought to a stop.
- Only after making sure that the rear legs are in contact with the bumper of the vehicle, engage the green lever to release the rear legs, continuing to push towards the inside of the vehicle compartment.

⚠ During these operations, a part of the weight of the assembly will rest on the operator, who must therefore be able to support and accompany all device movements.

- Lock the stretcher on the vehicle using the Spencer fixing system installed on the vehicle.

⚠ WARNING: DURING NORMAL STRETCHER MOVEMENTS, NEVER ENGAGE THE LEG LOCKING LEVERS, AS DOING SO COULD CAUSE THE STRETCHER TO FALL TO THE GROUND. THESE LEVERS MUST ONLY BE ENGAGED TO LOAD THE STRETCHER ON THE VEHICLE OR TO REACH THE INTERMEDIATE HEIGHT AS DESCRIBED ABOVE.

IT

EN

DE

FR

ES

PT

BG

HU

LT

EL

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

To unload the stretcher from the vehicle:

The intended fixing system is Z500900.

- Pull the stretcher towards the outside of the vehicle, gripping the rear part of the frame near the footrest. Support the weight of the stretcher until you feel the rear leg locking system engage. **As with the loading phase, the operator must be able to support the weight of the device.**
-  Do not grab any area not intended for handling as this may cause damage to the operator, the patient or the device.
- Pull the stretcher outwards until the front legs are fully open. **Do not lower the loading trolley from the loading surface until you have checked that the front leg locking system is properly engaged.**
- Finish extracting the stretcher from the vehicle.

11.7 ADJUSTING THE STRETCHER LENGTH

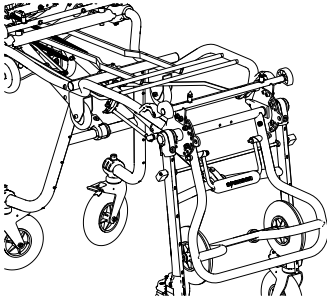
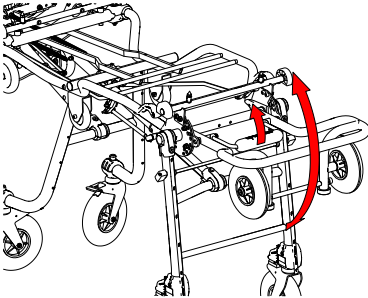
The maximum length of the stretcher can be reduced as follows to facilitate bringing the Spencer Cross stretcher in especially confined spaces (i.e. lifts):

- Hold the loading trolley securely and unlock it by pressing down or lifting the unlocking lever as shown in the adjacent image.
- Slowly lower the loading trolley with one hand, pushing it until it rests on the front legs.
- To return the trolley to the horizontal position, unlock it using the unlocking lever, supporting it as it raises and, again using the lever, align the trolley with the dedicated stops.



11.8 METHOD OF USE FOR COFFIN SUPPORT

If the intention is to use the stretcher to move coffins, place the positioning elements as indicated below.

FOOT SIDE SUPPORT

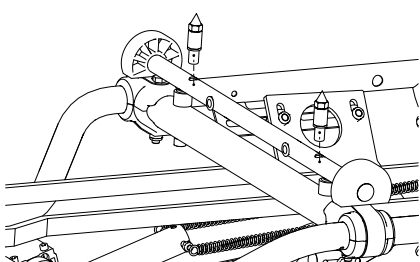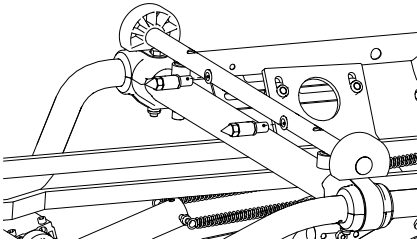HEAD SIDE SUPPORT

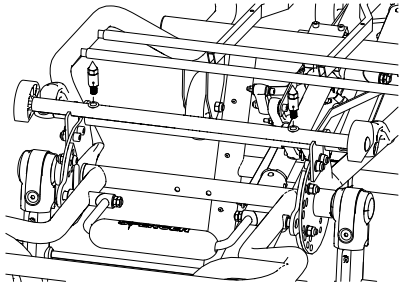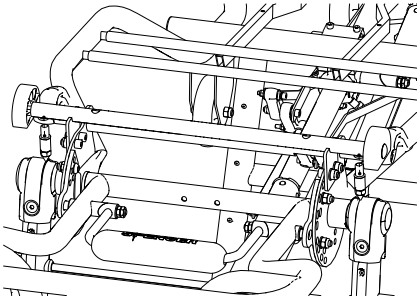12. CLEANING AND MAINTENANCE

Spencer Italia S.r.l. declines all responsibility for any direct or indirect damage which is the consequence of improper use of the product and spare parts and/or in any case of any repair carried out by a person other than the Manufacturer, who uses internal and external technicians authorised; to do so; moreover, doing so will invalidate the warranty.

- The operator must wear suitable personal protective equipment, such as gloves, goggles, etc. during all checking, maintenance and cleaning operations.
- Establish a maintenance schedule, periodic inspections and extend the average life span, if foreseen by the Manufacturer in the User Manual, identifying a reference person who meets the basic requirements set forth in the User Manual.
- The frequency of inspections is determined by factors such as legal requirements, type of use, frequency of use, and environmental conditions during use and storage.
- Repairs of products must be carried out by Spencer Italia S.r.l., who shall make use of specialised internal or external technicians who, using original spare parts, shall provide quality repair service in strict compliance with the technical specifications indicated by the Manufacturer. Spencer Italia S.r.l. declines any responsibility for any direct or indirect damage resulting from improper use of spare parts and/or any repair work carried out by unauthorised parties.
- Use only original or Spencer Italia S.r.l. approved components/spare parts and/or accessories in order to carry out all operations without causing alterations or modifications to the product.
- All maintenance and overhaul activities must be recorded and documented with the relevant technical operation reports. This documentation must be kept for at least 10 years after the end of the product's life and must be made available to the competent authorities and/or the Manufacturer when requested.
- Cleaning, provided for reusable products, must be carried out in accordance with the Manufacturer's instructions in the User Manual in order to avoid the risk of cross-infection due to the presence of body fluids and/or residues.
- The product and all its components must be washed and left to dry completely before storage.
- If the product requires lubrication, this must be done after cleaning and complete drying.

12.1 CLEANING

Failure to carry out the correct cleaning operations could increase the risk of cross-infection due to presence of body fluids and/or residues.
The operator must wear suitable personal protective equipment, such as gloves, goggles, etc. during all checking and cleaning operations.

Metal parts exposed to external agents undergo surface treatments and/or coating in order to obtain better resistance. Clean the exposed parts with water and delicate soap. **Never use solvents or stain removers. Do not use cleaning agents containing sodium hypochlorite as corrosion of components may occur.**
Rinse thoroughly with lukewarm water, making sure you have removed all traces of soap, which may deteriorate or compromise conditions and durability. **Avoid using high pressure water**, as it penetrates the joints and removes lubricant, creating the risk of corrosion on components. Let dry completely before storing. Drying after washing or after use in a wet environment must be natural and not forced. Do not use flames or other direct heat sources.
If **disinfecting**, use products that do not have a solvent or corrosive action on materials constituting the device. Be sure to take all precautions to ensure that there is no risk of cross-infection or contamination of operators.

12.2 ORDINARY MAINTENANCE

Establish a maintenance program and periodic checks, identifying a designated reference person. The person entrusted with device maintenance must ensure the basic requirements envisaged by the manufacturer in the following paragraphs.
All routine and special maintenance activities and all general overhauls must be recorded and documented with the relevant technical intervention reports. This documentation must be kept for at least 10 years after the end of the device's life and must be made available to the competent authorities and/or the Manufacturer when requested.

Routine device maintenance must be carried out by operators with specific qualifications, training and instruction in the use and maintenance of the device.

The operator must wear suitable personal protective equipment, such as gloves, goggles, etc. during all checking, maintenance and cleaning operations.

To ensure product traceability and protect the maintenance and servicing procedures for your devices, Spencer has made the following portal available to you www.spencer.it/en/services/report-an-issue

Checks to be carried out before and after each use and at deadline indicated above, are as follows:

- General functionality of the device
- Cleanliness of the device (remember that the failure of cleaning may cause the risk of cross infections)
- Correct tightening of all nuts, bolts and screws
- Absence of cuts, holes, tears or abrasions on the entire structure
- None of the tubes or metal sheets present bends or cracks
- Welded areas are intact, without any cracks or breaks
- Moving parts, wheels, levers, handles are intact and work properly
- Lubrication of moving parts
- Conditions of use of wheels and brake system
- The wheels are correctly fixed, are stable and turn properly
- The wheels are free of debris
- The device opens and locks properly
- The device opens and closes properly
- Functioning of springs
- The stretcher easily enters the vehicle
- The vehicle is equipped with a Spencer fixing system dedicated to the stretcher
- The coupling between the fixing system and the stretcher is suitable to guarantee the safety of the fixing

The frequency of inspections is determined by factors such as legal requirements, type of use, frequency of use, and environmental conditions during use and storage.

Please note that cleaning must be carried out as described in this manual and that functionality must be verified before and after each use. Spencer Italia S.r.l. declines any responsibility for improper functioning or damages caused to the operator/s by use of devices not subject to routine maintenance, invalidating the warranty.

Use only original or Spencer Italia S.r.l. approved components/replacement parts and/or accessories to carry out any operation without causing alterations or modifications to the device. Otherwise, we decline all responsibility regarding incorrect operation or any damage caused by the device to the operator, invalidating the warranty.

12.3 PERIODIC OVERHAUL

The device must be overhauled every year by the Manufacturer, who uses internal and external technicians specialised and authorized by the Manufacturer itself.

In the absence of an above-mentioned overhaul, the device must be **DECOMMISSIONED**, as the device no longer meets the safety requirements guaranteed by the Manufacturer at the time of delivery.

Spencer Italia srl declines any responsibility for incorrect operation or damage caused by use of devices that have not been regularly serviced.

The gas spring on the front legs must be replaced every 2 years. More heavy-duty or particular conditions of use may, however, require it to be replaced earlier. This need is evaluated by the manufacturer or by the service centre that carries out the overhaul. Failure to replace it within the prescribed period could compromise the safety of the stretcher.

Only overhaul activities carried out by specialised technicians authorised by the Manufacturer are considered valid by Spencer Italia S.r.l.

IT

EN

DE

FR

ES

PT

BG

HU

LT

EL

12.4 SPECIAL MAINTENANCE

Special maintenance can only be carried out by the Manufacturer, who uses internal and external technicians specialised and authorised by the Manufacturer itself.

Only maintenance activities carried out by specialised technicians authorised by the Manufacturer are considered valid by Spencer Italia S.r.l..
The end-user can replace only the spare parts indicated in § 15.

12.5 LIFE SPAN

The device, if used as described in the following instructions, has a life span of 5 years from the date of purchase, which can be extended following annual overhauls.

Overhauls must be carried out by the Manufacturer, who uses internal and external technicians specialised and authorised by the Manufacturer. In the absence of such annual overhauls, the device must be DISPOSED OF ACCORDING TO THE INFORMATION IN PARAGRAPH 16 AND NOTIFICATION MUST BE GIVEN TO THE MANUFACTURER.

The life span may be extended at the sole discretion of the Manufacturer or authorised centre if the safety requirements of the device are still met.

Spencer Italia S.r.l. will accept no responsibility for incorrect operation or damage caused by the use of devices that have not been overhauled by the Manufacturer or authorized centre, or that have exceeded the maximum allowable life span.

13. TROUBLESHOOTING TABLE

PROBLEM	CAUSE	REMEDY
The leg unlocking levers don't work properly or are difficult to operate.	The moving mechanisms are damaged.	Immediately remove the device from service and contact customer service.
	Means of connection between components have been lost.	
The stretcher is not correctly attached to the fixing system.	Wear or damage of the components that constitute the locking mechanisms.	Immediately remove the device from service and contact customer service.
	The fifth wheel has not been properly inserted in the fixing system.	Position the stretcher properly and check that the fifth wheel enters the dedicated part on the fixing system.
Damage to the structure	Improper use	Immediately remove the device from service and contact customer service.
	The moving mechanisms are damaged.	Immediately remove the device from service and contact customer service.
It is not possible to set the stretcher in intermediate position.	Something is interfering with the moving systems.	Make sure that nothing is interfering with the mechanisms.
	Levers have not been properly operated.	Carefully follow the instructions to place the stretcher in the intermediate position.
The front legs do not lock when unloading from the vehicle.	The moving mechanisms are damaged.	Immediately remove the device from service and contact customer service.
	The height of the loading floor is not suited to the device, the safety height is not being respected.	Adjust the loading floor height to meet the requirements described in this manual. If the loading floor does not allow adjustments, immediately remove the device from service and contact the service centre.
When the stretcher release lever is operated, the stretcher does not move from the fixing system and remains attached.	The magnet controlling the release, located inside the fifth wheel holder, has moved or been lost.	Manually unlock the stretcher and finish the unloading operation. At the end of service, check the positioning of the magnet, then restore its original position. If the problem persists, immediately remove the device from service and contact the service centre.
	The stretcher is used with a fastener that does not provide automatic release.	If necessary, request a different fixing system.

14. ACCESSORIES

ZS00900	HUX MORTUARY STRETCHER FIXING SYSTEM
CR90016	EXTENSION FOR MOVEMENT CROSS/STAR-K/AGILE
CR90024	LOWERING EXTENSION KIT CROSS/STARK/AGILE SERIES

15. SPARE PARTS

ZS00403	MOD.402 TWO-PIECES BLACK BELT BUCKLE IN DERLIN
---------	--

16. DISPOSAL

When devices and their accessories are no longer suitable for use, they can be disposed of as normal municipal solid waste if they have not been contaminated by special agents. Otherwise, follow the regulations in force regarding disposal.

Notice

The information contained in this document is subject to change without notice and is to be intended as a commitment by Spencer Italia S.r.l. subject to change.
The images are included as examples and may vary slightly from the actual device.

© Copyright Spencer Italia S.r.l.

All rights reserved. No part of the document may be photocopied, reproduced or translated into another language without prior written consent from Spencer Italia S.r.l.

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