



LIGHTWEIGHT STAIR CHAIR

4BELL SILVER+

MAIN FEATURES

- Frame in aluminum makes the chair lightweight and resistant at the same time
- Two telescopic handles at the front and two foldable ones behind the backrest help lifting and transport operations
- Rear wheels' locking mechanism makes sure the chair does not slide on a flat surface
- Semirigid seat and soft backrest are built to increase patient comfort
- Foldable footrest increases patient stability
- Easy to clean PVC seat and backrest are resistant to chemicals
- Pivoting front wheels makes the chair more maneuverable.



Class I Medical Device compliant with MDR (EU) 2017/745

— Tested according to —



PRODUCT DESCRIPTION: The 4BELL SILVER+ is a stair chair designed to transfer or evacuate patients with mobility limitations from multi-storey buildings. The chair is FDA compliant and built in accordance to the specifications set by EN 1865-1. If the stair chair is installed on the 4Bell Max+ bracket, it will comply with 10G crash test according EN 1789 standard. The frame of the chair is made in aluminum for increased resistance. Its total weight is 9,4 kg (20.7 lb).

INTENDED USE: Folding transport chairs like the 4BELL SILVER+ are intended for moving patients and positioning them in a seated position within the ambulance but are not designed for patient transport inside the ambulance itself. They should not be used for parking or stationary positioning.

397 lb
180 kg
LOAD CAPACITY

20.7 lb
9,4 kg
LOW-WEIGHT


FOLDABLE CHAIR

EASY TO STORE

TECHNICAL DATA

Code	ST10500
Length with handles retracted	76 cm / 29.9 in
Length with handles extracted	112 cm / 44 in
Width	49 cm / 19.8 in
Height	90,5 cm / 35.6 in
Dimensions when folded	49 x 90,5 x 29 cm / 19.8 x 35.6 x 10.5 in
Weight	9,4 kg / 20.7 lb
Load capacity	180 kg / 397 lb
Front wheels	Ø 10 cm / 3.93 in
Rear wheels	Ø 15 cm / 5.90 in
Materials	Al, Steel, PU, PVC
Seat sheet color	

Note: Dimensions and weight are subject to tolerances of ± 1 cm (± 0.39 in) and ± 0.2 kg (± 0.44 lb) respectively.

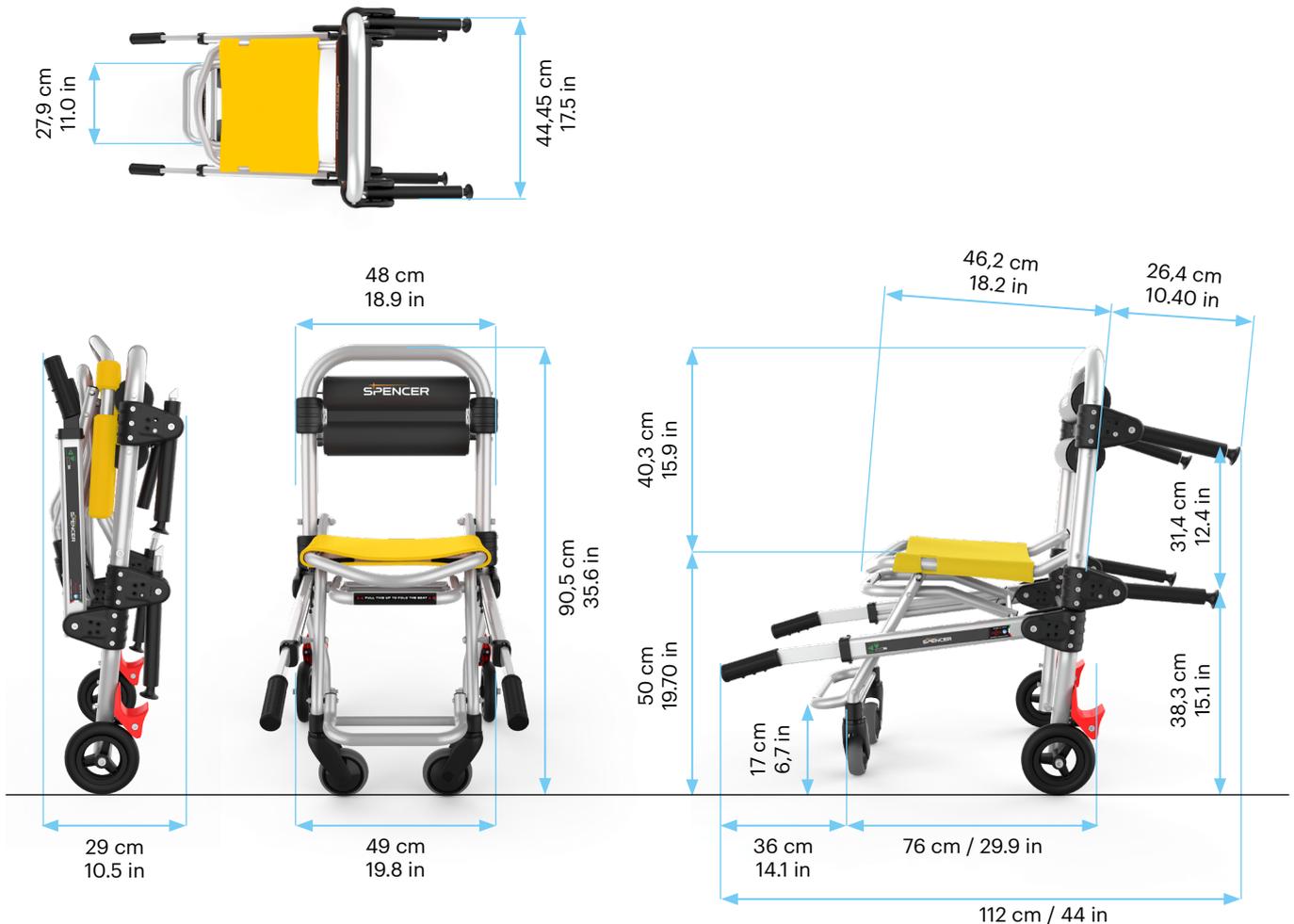



THE PRODUCT INCLUDES

- Restraint belts.

REQUIRED FASTENING SYSTEM



ST10453

4BELL+ MAX
10G tested fastening system

- It is installed on ambulances, special vehicles or aeroplanes to secure all 4Bell+ series chairs in the folded position.
- It can be fixed to either the wall or floor, depending on the structure and configuration of the vehicle.

ACCESSORIES



ST10452	Seat XL
ST10450	Armrests
ST10451	Rear auxiliary handles 4BELL+

Rev.0 (09/03/2026)
UNCERTIFIED COPY

©Copyright - This document is property of SPENCER ITALIA S.r.l and can't be reproduced, used and disclosed without the written approval of SPENCER ITALIA SRL.