



COMPACT STAIR CHAIR WITH TRACKS

4BELL STAIR+


MAIN FEATURES

- Lightweight and durable aluminum frame.
- Reinforced tracks ensure safe, controlled descent on stairs.
- Extendable rear handle offers enhanced control for stair descent.
- Telescopic handles aid in lifting over obstacles.
- Rear wheel locking mechanism prevents sliding on flat surfaces.
- Semi-rigid seat and soft backrest enhance patient comfort.
- Foldable footrest adds stability for the patient.
- Pivoting front wheels increase maneuverability.
- Easy-to-clean PVC seat and backrest resist chemical exposure.



Class I Medical Device compliant with MDR (EU) 2017/745

— Tested according to —

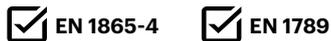


PRODUCT DESCRIPTION: The 4BELL STAIR+ is a foldable EMS transport chair specifically designed for first responders to safely transfer or evacuate patients with limited mobility from multi-story buildings. Equipped with dual tracks, the chair allows for controlled descent down stairs, minimizing the need for lifting. Constructed with a durable aluminum frame, it weighs just 14.7 kg (32.4 lb), offering both strength and ease of handling. It is 10G tested when fastened to its dedicated anchoring system.

INTENDED USE: Folding transport chairs like the 4BELL STAIR+ are intended for moving patients and positioning them in a seated position within the ambulance but are not designed for patient transport inside the ambulance itself. They should not be used for parking or stationary positioning.

397 lb
180 kg

LOAD CAPACITY

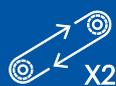

STAIR TRACKS



MULTI-HEIGHT EXTENSION HANDLE



FOLDABLE CHAIR

TECHNICAL DATA

Code	ST10501
Length with handles retracted	76 cm / 29,9 in
Length with handles extracted	112 cm / 44,09 in
Length with tracks and telescopic handles extracted	128 cm / 50,3 in
Width	49 cm / 19,2 in
Height	94 cm / 37 in
Max. height with main handle fully extended	147,5 cm / 58.07 in
Dimensions when folded	49 x 94 x 29 cm / 19,2 x 37 x 11,4 in
Weight	14,7 kg / 32.4 lb
Load capacity	180 kg / 397 lb
Front wheels	Ø 10 cm / 3,93 in
Rear wheels	Ø 15 cm / 5,90 in
Materials	Al, Steel, PU, PVC
Seat sheet color	

Note: Dimensions and weight are subject to tolerances of ± 1 cm (± 0.39 in) and $\pm 0,2$ kg (± 0.44 lb) respectively.

THE PRODUCT INCLUDES

- Two-piece restraint belts.

ACCESSORIES



ST10455

4BELL STAIR+ MAX

10G tested fastening system



ST10452	Seat XL
ST10450	Armrests
ST10451	Rear auxiliary handles 4BELL+
ST10454	Headrest 4BELL STAIR+



Rev.0 (09/03/2026)
UNCERTIFIED COPY

©Copyright - This document is property of SPENCER ITALIA S.r.l and can't be reproduced, used and disclosed without the written approval of SPENCER ITALIA SRL.

Spencer Italia s.r.l. a socio unico soggetta a direzione e coordinamento da Protect Medical Holding GmbH
Via Provinciale, n° 12 - 43038 Sala Baganza (PR) Italy - +39 0521 541111 - info@spencer.it - www.spencer.it