

CARRYING SHEET WITH 10 HANDLES AND INTEGRATED RESTRAINTS BELTS

# PAX WOW

## MAIN FEATURES

- Carrying sheet with increased load capacity for transporting patients up to 250 kg.
- Removable aluminium rods to increase sheet stiffness allowing load distribution over a larger area and limited flexing.
- Made of highly durable and non-absorbent materials, it ensures maximum safety and hygiene during use.
- Two adjustable shoulder straps to improve load distribution and reduce operator effort.
- Designed to ensure operation even with only two operators\*
- Patient handles system for improved grip and increased sense of security.
- Footrest adjustable in two positions to enhance longitudinal containment and prevent slipping.
- Integrated belts to ensure to ensure patient safety.
- Easy to fold, occupies minimal space when rolled up.



*\*The number of operators must be proportional to the weight to be lifted, in accordance with occupational safety regulations.*



MANUFACTURED BY X-CEN-TEK



Compliant with EN 1865-1 standard

Class I medical device compliant with Regulation (EU) 2017/745 (MDR)

**PAX WOW** is a carrying sheet designed for the recovery and transport of patients with no signs of spinal trauma, from the scene of first response to a stretcher or rescue vehicle. It is particularly suited for stairway descent, where rapid intervention or confined spaces prevent the use of stretchers or evacuation chairs.

Its special anatomic shape, wider in the head-shoulder area provides enhanced lateral containment, improving patient stability and safety during transfer and transport. It features a three-segment foldable aluminium support system, connected by elastic polyester braid, which can be inserted into dedicated lateral guides along the sides of the sheet. The use of these lateral supports, which are self-assembling and easy to insert, allows the full length of the device to be used, providing reliable support for the patient.

Equipped with a top pocket with magnetic closure, it allows the insertion of a pillow or blanket for patient comfort and, after use, storage of rods and shoulder strap, keeping accessories always within reach and minimising bulk. The sheet design enables transport of the patient in both supine and seated positions, adapting to different operational needs.



**TECHNICAL DATA**

<b>Code</b>	<b>284970201</b>
Length	256 cm
Width	Max 96 cm (top), min 70 cm (footside)
Wrapped storage (1)	47 x 30 x 14 cm
Shoulder strap adjustment positions	6
Reinforced handles	10
Shoulder strap width	5 cm
Length of upper/lower handles	15 / 26 cm
Rod length	115 cm
Dimensions of single folded rod	7,4 x 41 x 3 cm
Weight	approx. 5 kg
Load capacity	250 kg
Materials	Nylon, PVC, Aluminium
Color	


(1) The measurement refers to the sheet when rolled up. It is purely indicative, as the device can be folded or rolled in different ways. Dimensions and weight have a tolerance of  $\pm 1$  cm and  $\pm 0,2$  kg respectively.



The 10 handles, made from reinforced nylon webbing, provide a secure grip, high strength, and long-lasting durability, ensuring safe and comfortable patient transport.



The footrest, adjustable in two positions, enhances patient containment, providing greater stability during transport.

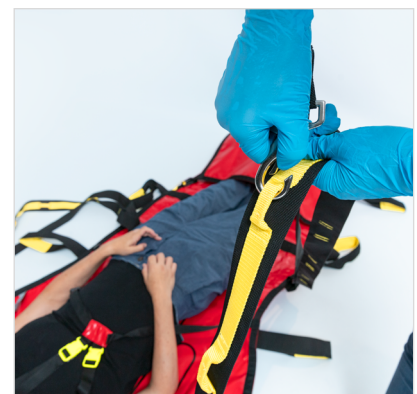

The adjustable lateral shoulder straps guide the rescuer to perform proper lifting movements, helping to significantly reduce muscle strain.



The aluminium rods increase the rigidity of the sheet, allowing the load to be distributed over a larger surface and limiting flexing.



The aluminium rods can be removed through dedicated convenient holes, folded into three sections, and stored in the built-in pocket within the sheet.



Integrated belts secure the patient safely, ensuring effective immobilisation during transport.

Dimensions of single folded rod: 7,4 x 41 x 3 cm



The inner handles are designed to provide the patient with enhanced security during transport.



Reinforced, foldable structure for minimal storage, combining compactness and ease of transport.



Rev.0 (09/02/2026) – UNCERTIFIED COPY  
DISTRIBUTED BY SPENCER ITALIA S.R.L.

©Copyright - This document is property of SPENCER ITALIA S.r.l and can't be reproduced, used and disclosed without the written approval of SPENCER ITALIA SRL.