




SELF-LOADING INCUBATOR TROLLEY

# CROSS INCUBATOR

## MAIN FEATURES

- Stainless steel stretcher with micro shot-peening for enhanced longevity
- Shortenable frame for improved mobility in confined spaces
- The brakes on the rear wheels ensure greater stability and control during parking.
- Colour-coded levers make trolley loading into the ambulance more intuitive.
- Adjustable loading height system allows setting the correct height when changing vehicles or configurations.

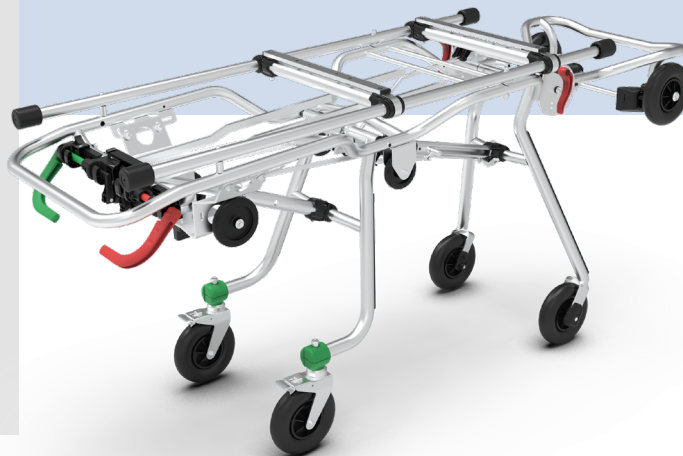

Tested according to EN 1789

Class I Medical Device compliant with MDR (EU) 2017/745

**PRODUCT DESCRIPTION:** The CROSS INCUBATOR is a stainless steel, self-loading trolley with adjustable heights for the efficient and safe transport of incubators. The legs fold independently and, by activating the dedicated controls during the unloading procedure, open automatically. The front trolley can be lowered or rotated upwards for better operation in confined spaces or small lifts.

330.7 lb  
150 kg  
LOAD  
CAPACITY

FOLDABLE  
FRAME

ADJUSTABLE  
LEVELS

## TECHNICAL DATA

Code	CR60010-ZO
Length	77.5 cm / 77.95 in
Length folded	166 cm / 65.35 in
Width	56 cm / 22 in
Weight	40,5 kg / 89.3 lb
Loading capacity	150 kg / 330.7 lb
Wheels	Ø20 cm / Ø7.87 in
Loading height <sup>(1)</sup>	From 58 to 67 cm / From 22.8 to 26.4 in
Frame material	Stainless steel

Dimensions and weights are subject to tolerances of  $\pm 1$  cm (0.39 in) and  $\pm 2$  kg (4.41 lb) respectively.  
 (1) Height measured from the ground to the loading wheels..

## ACCESSORIES

- Code EN90003 – END-T Instrument Table 10G for Stretcher
- Code ST42720 – S-MAX PRO 10G FIXATION SYSTEM

Rev.0 (10/06/2025)  
 UNCERTIFIED COPY

©Copyright - This document is property of SPENCER ITALIA S.r.l and can't be reproduced, used and disclosed without the written approval of SPENCER ITALIA SRL.