




SELF-LOADING INCUBATOR TROLLEY

# CROSSOVER INCUBATOR

## MAIN FEATURES

- Stainless steel stretcher with micro shot-peening for enhanced longevity
- The integrated Twist system allows the front wheels to be unlocked, making them swivel for improved manoeuvrability
- Shortenable frame for improved mobility in confined spaces
- The brakes on the rear wheels ensure greater stability and control during parking.
- Colour-coded levers make trolley loading into the ambulance more intuitive.
- Adjustable loading height system allows setting the correct height when changing vehicles or configurations.

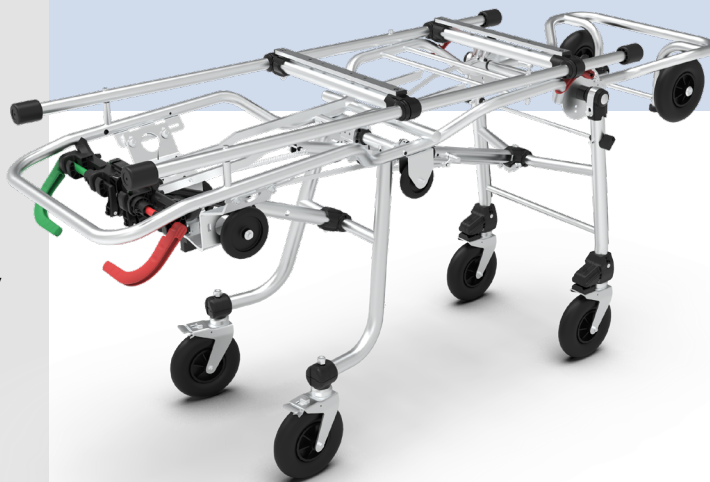

Tested according to EN 1789

Class I Medical Device compliant with MDR (EU) 2017/745

**PRODUCT DESCRIPTION:** The CROSSOVER INCUBATOR is a stainless steel, self-loading trolley with adjustable heights for the efficient and safe transport of incubators. The legs fold independently and, by activating the dedicated controls during the unloading procedure, open automatically. The front trolley can be lowered or rotated upwards for better operation in confined spaces or small lifts. The TWIST system, integrated into the trolley, unlocks the front wheels, allowing them to pivot for better maneuverability.

330.7 lb  
150 kg  
LOAD  
CAPACITY

FOLDABLE  
FRAME

3  
ADJUSTABLE  
LEVELS

## TECHNICAL DATA

Codes	CR60008-ZO	CR60009-ZO
Length	197 cm / 77.5 in	197 cm / 77.5 in
Length folded	166 cm / 65.35 in	166 cm / 65.35 in
Width	56 cm / 22 in	56 cm / 22 in
Weight	43,5 kg / 96 lb	43,5 kg / 96 lb
Loading capacity	150 kg / 330.7 lb	150 kg / 330.7 lb
Wheels	Ø20 cm / Ø7.87 in	Ø20 cm / Ø7.87 in
Loading height <sup>(1)</sup>	From 58 to 67 cm From 22.8 to 26.4 in	From 47 to 55 cm From 18.50 to 21.65 in
Frame material	Stainless steel	Stainless steel

*Dimensions and weights are subject to tolerances of ± 1 cm (0.39 in) and ± 2 kg (4.41 lb) respectively.  
 (1) Height measured from the ground to the loading wheels.*

## ACCESSORIES

- Code EN90003 – END-T Instrument Table 10G for Stretcher
- Code ST42720 S-MAX PRO 10G FIXATION SYSTEM
- Code ST0033 INCUBATOR BLOCKING SYSTEM DRÄGER TI500